RESMED

Important: Remove this guide before handing the device to the patient. Pressures should only be adjusted by qualified personnel.

S8 Compact™

QUICK SET-UP GUIDE

Set-up is quick and easy . . . Just Go!

The following settings have been pre-set to accommodate typical patient treatment values.

S8 Compact™

Altitude:	0–2000 ft (0–609 m)
Ramp:	20 minutes
Max Ramp:	45 minutes
Start CPAP:	4 cm H ₂ O
CPAP:	9 cm H ₂ O
Reminders:	OFF

1

Accessing the clinical menu and changing settings and options

- Turn on the device and wait until the standby (Ramp) screen appears on the LCD. At this stage you can adjust the ramp time.
- Press the **Down** and **Right** keys ( + ) simultaneously for 3 seconds to access the clinical menu.
- Use the **Up**  and **Down**  keys to scroll through the menu.
- Use the **Left** key  to change each setting as required.
- Press the **Right** (exit) key  to return to the Patient Menu screen.

2

Accessing and setting the reminder menu


- Turn on the device and wait until the standby (Ramp) screen appears on the LCD. At this stage you can adjust the ramp time.
- Press the **Left** and **Down** keys ( + ) simultaneously for 3 seconds to access the reminder menu.
- Use the **Up**  and **Down**  keys to scroll through the menu.
- Use the **Left** key  to change each setting as required.
- Press the **Right** (exit) key  to return to the Patient Menu screen.

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S8 Compact™

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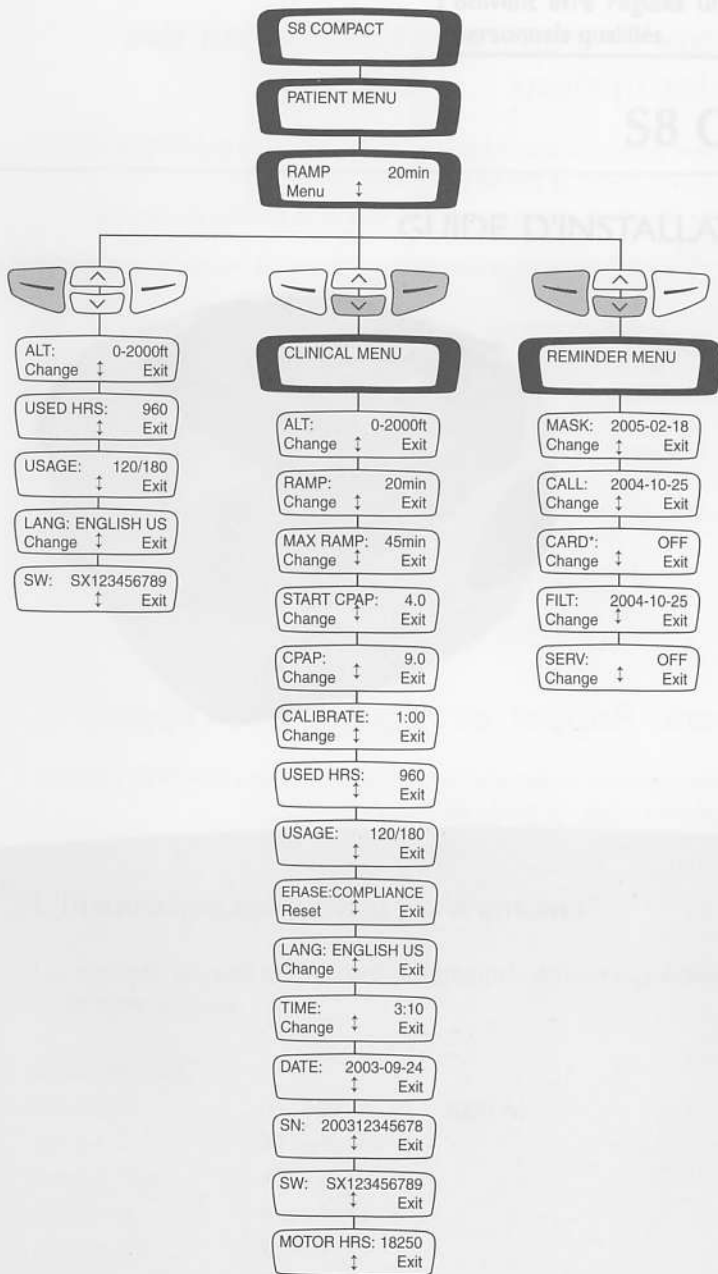* Feature not available on this device.